

Glycol Buffer - Instructions

LogTag Recorders' glycol buffer is used where vaccines need to be monitored during storage or transport. It mirrors the temperature environment inside a vaccine vial and is compliant with the CDC requirements for vaccine storage monitoring.

The buffer is used with LogTag®'s ST100K-15 or ST100K-30 sensors connected to any TRED30-7, PTED30-7 or TREX-8 temperature recorder and approx. 57ml (1.9 fl oz) of glycol. These items are not part of the kit and need to be ordered separately.

The glycol buffer contains three parts: the vial, the seal and the cap. To use the buffer with one of the sensors, please fill and seal the vial according to these instructions:



Glycol buffer (seal, cap and vial) with ST100K sensor



1. Fill the vial with glycol to the top mark.




2. Press the seal into the vial until the underside of the lip is flush with the top of the vial.




3. Screw on the cap until it is just touching the seal.
Do not tighten yet.



4. Insert the sensor into the cap's opening and slowly push until the blue heat shrink touches the cap. The sensor now reaches into the vial level with the bottom mark.

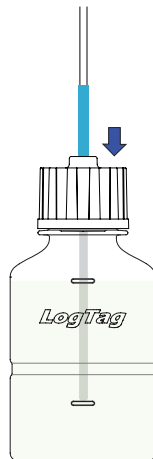

5. Firmly tighten the bottle cap without overtightening.

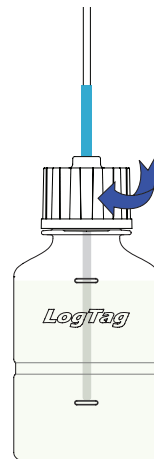

Connect the sensor to a compatible LogTag® recorder and place the glycol buffer with the vaccines to be monitored. If you follow these instructions, the sensor will stay in the center of the vial and the seal will not leak even when the buffer is tipped over.

Follow the instructions for configuring and using the recorders in the LogTag® Analyzer User Guide.



WessexPower

Solutions for High-tech Instrumentation

www.wessexpower.co.uk

 **LSTechnology**[®]

 E-commerce Instrumentation

www.loggershop.co.uk

Telephone: +44 (0) 1929 459 459 | Email: Sales@WPLS.co.uk